

LIPOZERO

Special Body & Face Line care



BSG AESTHETICS SOLUTIONS CO.

 EunSung

Contents


- 1) **What is Lipozero**
- 2) **Technology & Mode of Action**
 - **Multipolar RF**
 - **Vacuum Suction**
 - **Cavitation**
 - **Low Level Laser**
 - **LED**
- 3) **Component Overview**
 - **Main body**
 - **Left Handpiece**
 - **Right Handpiece**
- 4) **Clinical Indication**
- 5) **Summary**
- 6) **Specification**

1. What is Lipozero



LIPOZERO is the solution for total body-line management !

Two types of handpieces can treat from the face and body at the same time !

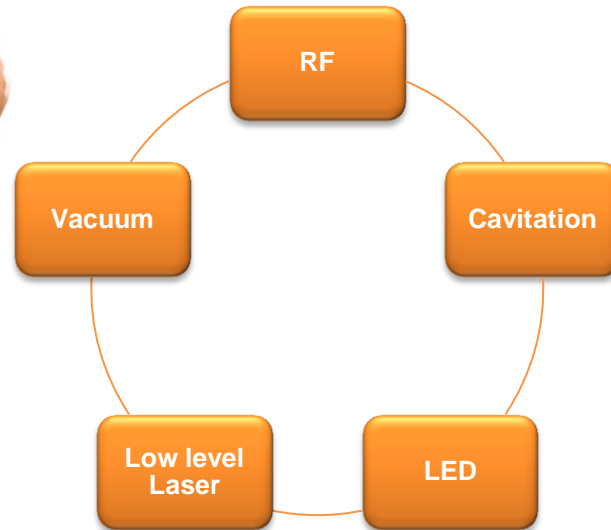

The combination of a different technologies for improved fat reduction.

Reduces the fat and cellulite at the same time, leading to improvement of over all healthier body.

Solves out the body care & Anti-aging at the same time.


***Professional Body & Face Rejuvenation
system in your clinic.***

2. Technology & Mode of Action

2.1. Multipolar RF

- Lipozero's Multipolar RF is one of the latest technological breakthrough in the fields of radiofrequency. Innovative technology harnesses the inherent benefits of RF yet overcomes the limitations of previous mono & bi-polar technologies.
- Completely safe and extremely efficacious. Optimal RF energy to the skin surface and subcutaneous fat layers resulting in immediate visible body and facial contouring effects without surgery or downtime.

- 1) Immediate collagen contracture & neoformation afterwards.
- 2) Increased fat cell metabolism.
- 3) Improvement of microcirculation & lymphatic drainage.

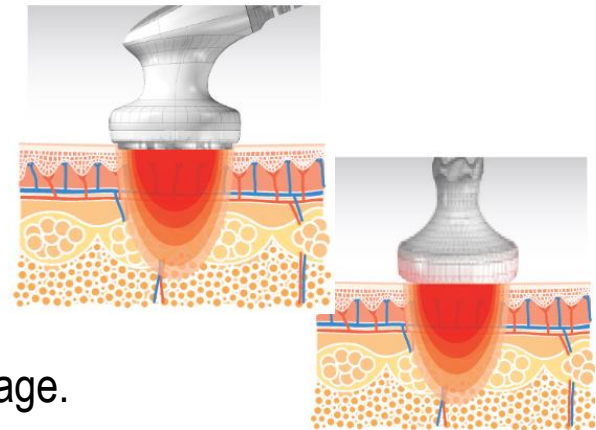


2. Technology & Mode of Action

2.2. Vacuum Suction

The suction effect occurs due to a reduced atmospheric pressure inside the applicator hole. As the air is sucked out, the skin and subcutaneous tissue is lifted upwards, folds then compresses.



- Increased local blood circulation & stimulates lymphatic drainage.
- Facilitates fibroblast activity & promotes vasodilation and diffuses oxygen.
- Enables precise & more up-close delivery of energy.

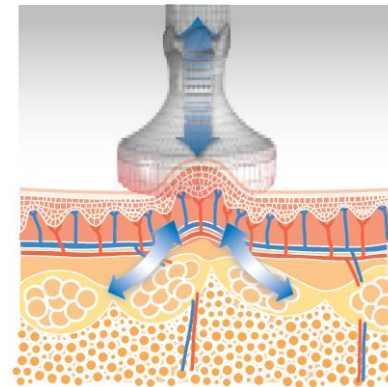

2. Technology & Mode of Action

2.3. Cavitation


- 34KHz uniquely engineered flat ultrasound technology.
- Acoustic thermal & acoustic mechanical effects heat & vibrate the treatment area and stimulate the fat cells to release their lipid content.

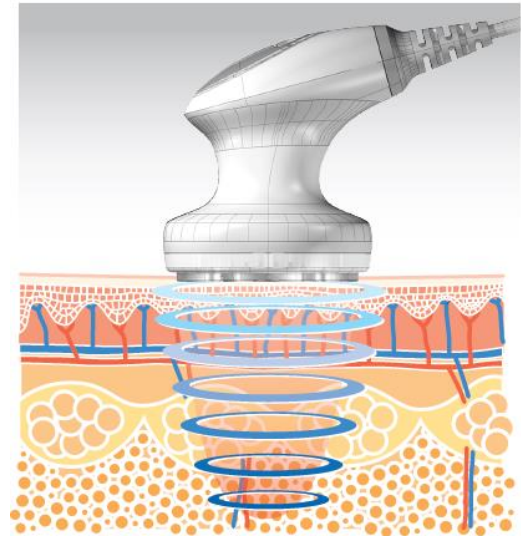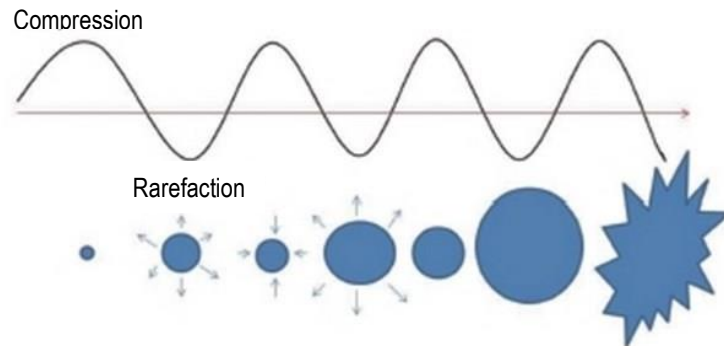

2. Technology & Mode of Action

2.4. Low Level laser



- 650nm low level laser.
- Photochemical biostimulation is generated by targeting the chromophore of Cytochrome C Oxidase.
 - 1) Disrupt the fat cell by creating transient pore at adipocyte cell membrane.
 - 2) Break down the triglycerides to 3 free fatty acid & glycerol by stimulating the lipase.

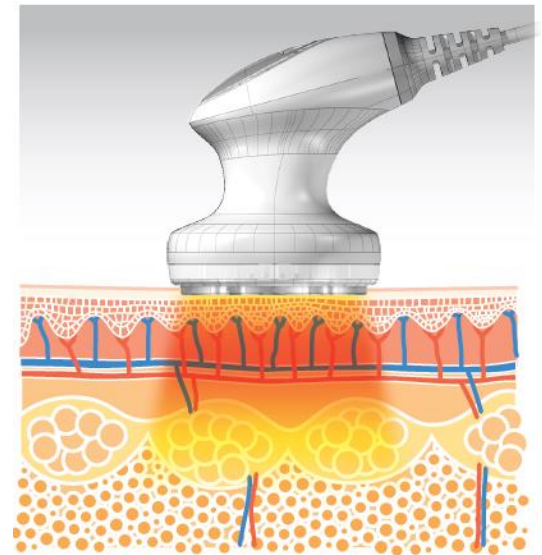

2. Technology & Mode of Action

2.5. LED

- Blue(415nm) & Red(620nm) LED therapy.
 - 1) Blue(415nm) works by killing the bacteria deep in the skin that is responsible for acne.
 - 2) Red(620nm) stimulates the fibroblasts that produce collagen, which gives young skin its plump look.

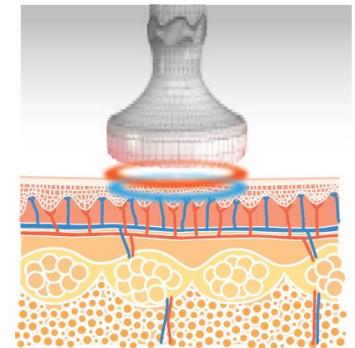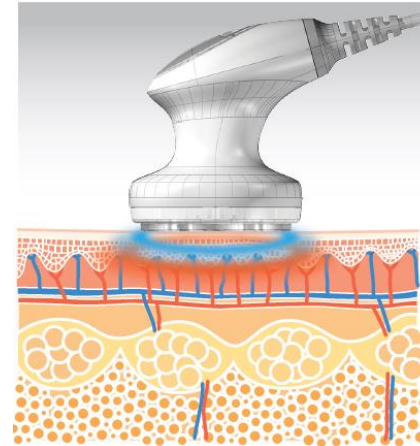

3. Component Overview

3.1. Main Body




3. Component Overview

3.2. Left Handpiece

Four in One easy to use handpiece. RF / Ultrasound / Low Level Laser / LED




[RF]

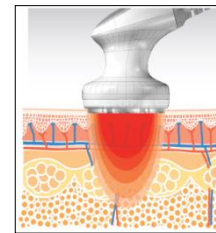


0.8 MHz

[Cavitation]

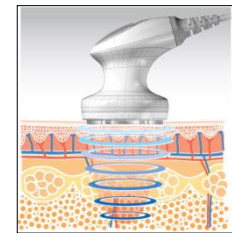


35 – 45 kHz

[Low level Laser]

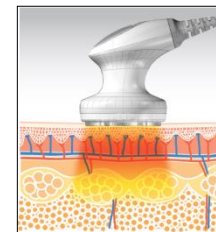

650 – 660 nm

[LED]

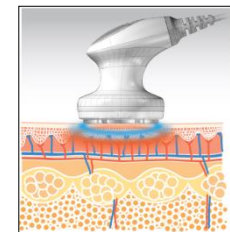


- Blue (415 nm)
- Red (620 nm)
- Blue & Red


3. Component Overview

3.3. Right Handpiece

Three in One easy to use handpiece. RF / Vacuum Suction / LED

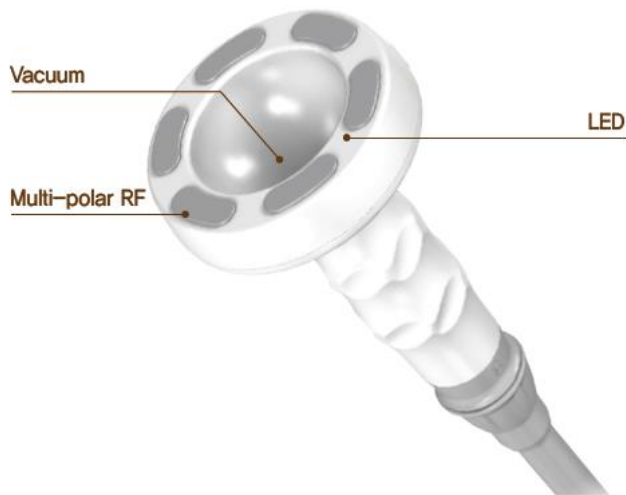


[RF]

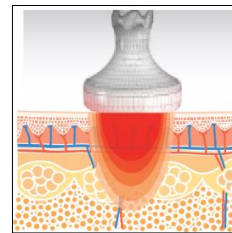


0.8 MHz

[Vacuum]

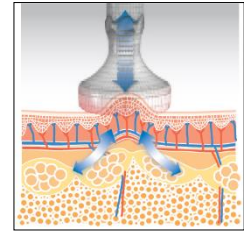

4 Vibrating Pulses

[LED]

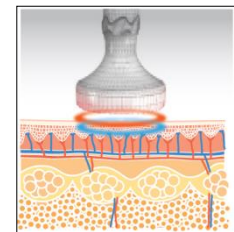

- Red (620 nm)
- Large & Medium
- Blue (415 nm)
- Small

For Body



Large

For Face & Body



Medium

For Face



Small

4. Clinical Indication

- Improvement of body-shape
- Cellulite improvement
- Post-liposuction
- Skin rejuvenation
- Improvement of facial-line
- Volumetric dermal therapy
- Wrinkle treatment

5. Summary

- Treatment Versatility. Multiple applicators and program settings allow you to tailor each treatment to the needs of your customers - from more aggressive therapy to a gentle postoperative session after liposuction.
- 5 Proven body contouring & skin care technologies - **RF, Cavitation, Vacuum, Low Level Laser, LED** - are in one device.
- Multifunctional device yet easy to operate.
- No consumables or disposables.
- Pain-Free & No Downtime.
- Engineered with proven body-contouring technology. Safe & Efficacious treatment guaranteed.

6. Specification



Technologies	RF, Vacuum Suction, Cavitation, Low Level Laser, LED
Power	AC 100 – 240 V, 50/60 Hz
Radiofrequency	0.8 MHz
Cavitation	35 – 45 kHz
Low Level Laser	650 – 660nm
LED	Blue (415nm), Red (620nm)
Dimension (W x D x H)	360 x 415 x 225 (mm)
Weight	8.5 kg



BSG AESTHETICS SOLUTIONS CO.

Thank you!

FOR FURTHER INFORMATION, PLEASE CONTACT US
AT 213.536.7229
or email us at
info@biospecgrp.com